



Ophthalmic Microscope

Hans Heiss System

HHSM-6622

Hans Heiss



MERCOFRAMES
OPTICAL CORP



info@mercoframes.com
alan@mercoframes.com



www.mercoframes.com



5555 NW 74 AVE.
MIAMI, FL. 33166



P: 305-882-0120
F: 305-882-0140

Features:

- Ali durable metal materials.
- Anti-mildew and anti-moisture optical lens.
- Simple installation and operation, no focus adjustment.
- Compact and generous design.
- Light weight for easy carrying.
- Workable with Stereo Indirect Ophthalmoscope to meet different posterior surgery requirements.

Image Inverter

Image Inverter is a kind of Adapter, used to put between Objective Lens and Eyepieces Lens of an Operating Microscope to provide an upright image for observer, while he is doing retina/vitreous surgery.

The image can be inverted by angle of 180°.

There is a thick metal stick to open and close inverting by moving up and down, without any focus adjustment. The operation is quite easy and convenient, even for beginners.


It is applicable in ZEISS, TOPCON, MOLLER, ALCON, TAKAGI, INAMI, LEICA & VISUMATE etc. Ophthalmic Operating Microscopes.



Front




Side

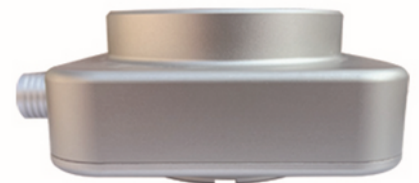

Back

Lens Set in details:

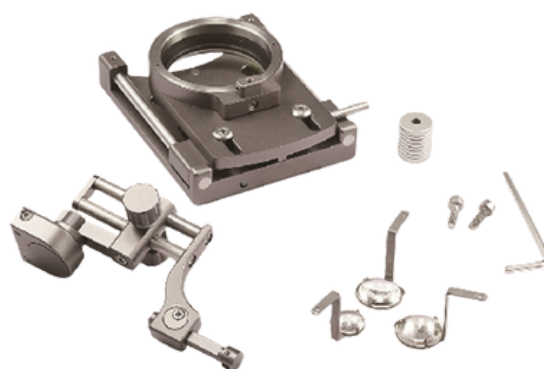

Technical Specifications:

Lens (3 pieces)	Field of View (F175mm/F200mm)	Image Magnification	Lens Diameter
No.1 Lens	102°/120°	0.34X	19mm
No.2 Lens	95°/112°	0.42X	13mm
No.3 Lens	80°/95°	0.74X	19mm

Related Matching Models

Brands	Models
Leica	M500, M501, M520, M620, M650, M651, M655, M690, M695, M820, M822, M840, M841, M844.
Zeiss	Lumera, Lumera T, MDI, MDO, MDU, OPMII, OPMI6, OPMI6-CFC, OPMI6S, OPMI6-SFC, OPMI6-SFR, OPMI6-CH, OPMIPico, OPMIVa- no, OPMIVISU 150, OPMIVISU 160, PMIVISU 200, OPMIVISU 210.
Topcon	OM 800 series, OMS 710, OMS 90.
Moller-wedel	Hi-R 900
Takagi	OM-18, OM-9
Alcon	LuxOR (Illumin-i)
Inami	L-0980A
Prescott	OMNI- FLEX
Huvitz	HOM-700
Visu Mate	SM-2000L, SM-3000L



MERCOFRAMES
OPTICAL CORP



info@mercoframes.com



alan@mercoframes.com



www.mercoframes.com

Hans Heiss



5555 NW 74 AVE.
MIAMI. FL. 33166



P: 305-882-0120
F: 305-882-0140